



SGX-E20R

HANDPIECE FOR ZYGOMATIC IMPLANTS

Easy and efficient solution for Zygomatic Implant

SGX-E20R can make your Zygomatic implant treatment easier and more efficient. Benefit from the unique shape specifically designed for Zygomatic implant treatment. NSK's special sealing system at the tip of the instrument prevents contamination from entering the instrument and the smooth surface makes decontamination easy.

The Optimal angle for Zygomatic implants


Compared to a conventional surgical contra-angle, SGX-E20R is specifically angulated to optimise access to the Zygomatic bone and greatly improves visibility.





Efficient cooling by external water supply

The external irrigation nozzle of the SGX-E20R cools the treatment field and bur efficiently.

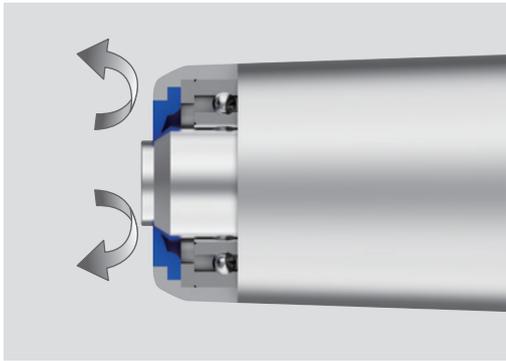

Original sealing system

Preventing contaminants from entering the instrument makes decontamination easier and prolongs the life of the instrument.



Gear Ratio	MODEL	ORDER CODE
20:1 Reduction	SGX-E20R	HA1200

- Stainless Steel Body • Max. Speed : 2,000 min⁻¹ • Max. Torque : 50 Ncm
- Twist Chuck • Angle Handpiece • External cooling
- For Osteotomy preparation / insertion of Zygoma implants for Maxillofacial / oral surgery



Washable in the thermodisinfectant



Autoclavable